
LYKKELE


Jilt

IPL Hair Removal Device

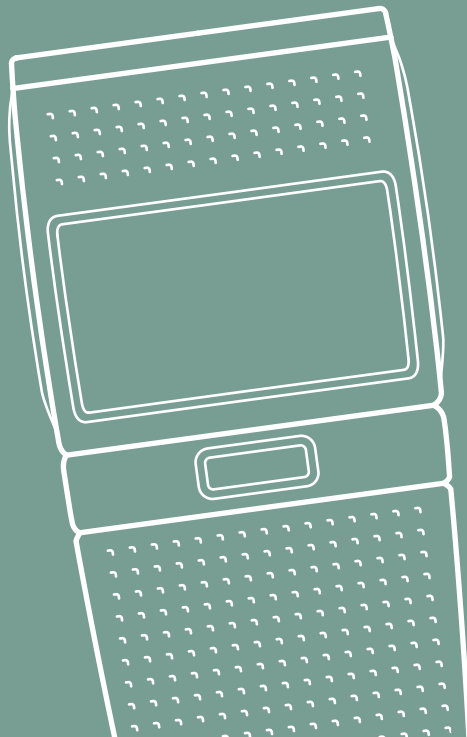Jilt

IPL Hair Removal Device

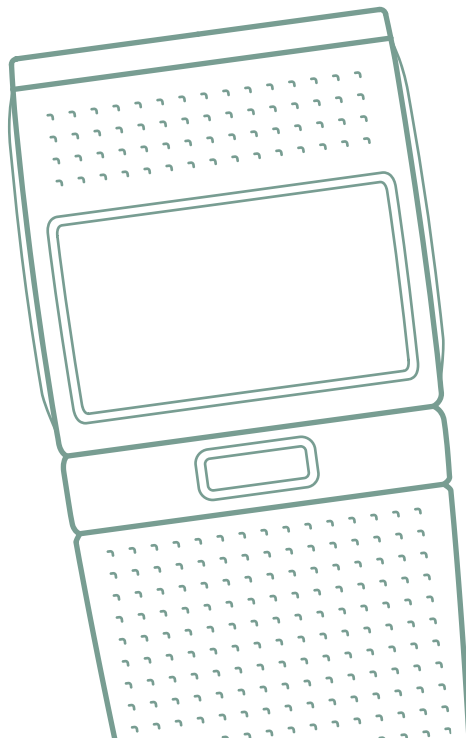Jilt

IPL Hair Removal Device

Before you use the IPL hair removal device, read the manual. Do not throw the manual out.

This is the manual's short version. To find the manual's full version in different languages, scan this QR code or visit: <https://www.lykkele.pl/do-pobrania>

Safety information

- Use the IPL hair removal device only as we indicate in this manual.
- Do not turn the IPL hair removal device on if it's damaged – e.g., it doesn't work properly, its case or its power adapter's plug is broken.
- If the IPL hair removal device is damaged, send it to us for repair. Otherwise, your warranty expires. We are not responsible for any damage that an IPL hair removal device which wasn't repaired in our service may cause.
- Keep the IPL hair removal device out of children's reach.
- Do not use the IPL hair removal device near water.
- If the IPL hair removal device falls into the water, do not touch it. Unplug its power adapter immediately.
- Do not put the IPL hair removal device under water or any other liquid.
- Do not touch the IPL hair removal device with wet hands.
- The IPL hair removal device emits light – don't look at it directly.
- Do not use the IPL hair removal device anywhere near your eyes.
- If you're not using the IPL hair removal device – turn it off and unplug its power adapter.
- Before you clean the IPL hair removal device or its power adapter, unplug it. Clean them with a slightly damp cloth.
- When you unplug the power adapter, hold onto the plug, not the cord. Do not touch the plug with wet hands.

Who shouldn't use the IPL hair removal device:

- people younger than 18;
- men on their facial hair; and
- people who have dark skin. If you use the IPL hair removal device on dark skin, it may cause burns, blisters, or skin discoloration.

Do not use the IPL hair removal device if:

- you have freshly tanned skin; or
- your skin was recently exposed to sunlight.

If you use the IPL hair removal device on such skin, it may cause burns, blisters, or skin discoloration.

Safety information

Do not use the IPL hair removal device if:

- you're pregnant or breastfeeding;
- you've had a skin grafting surgery less than four weeks ago;
- you have eczema, psoriasis, skin injuries, open wounds, or skin infections such as cold sores. In this case, you may use the IPL hair removal device after your skin lesions have healed;
- you have diabetes-related skin problems;
- you have epilepsy;
- you have a pacemaker, an incontinence device, or you use an insulin pump;
- you have gonorrhea, polymorphic light eruption, urticaria, lupus erythematosus, etc.;
- you have skin cancer or a doctor suspects that you might have it;
- you've had radiation therapy or chemotherapy in the last 3 months;
- you've taken isotretinoin in the last 6 months;
- you're currently using:
 - fruit acids,
 - beta-hydroxy acids,
 - all-trans-retinoic acid,
 - retinoids,
 - azelaic acid,
 - tetracycline,
 - phenothiazine,
 - thiazide,
 - diuretics,
 - sulfonylureas,
 - sulfonamides,
 - trimethoprim,
 - fluorouracil,
 - vinblastine,
 - griseofulvin,
 - painkillers or anti-inflammatory drugs, such as: Aspirin, Ibuprofen, or Acetaminophen.

Tip | People with disabilities should use the IPL hair removal device with a guardian present if they can't use it on their own.

Tip | If you're not sure if you may use the IPL hair removal device, contact your family doctor, surgeon, or dermatologist.

What you'll find in the box

- IPL Hair Removal Device
- Power adapter
- Shaving razor
- Safety glasses
- User manual

Before you use the IPL hair removal device

Tip | Do not expose your skin to sunlight or ultraviolet radiation for 4 weeks before you use the IPL hair removal device.

Skin						
Hair						
	✓	✓	✓	✓	✗	✗
	✓	✓	✓	✓	✗	✗
	✓	✓	✓	✓	✗	✗
	✓	✓	✓	✓	✗	✗
	✗	✗	✗	✗	✗	✗
	✗	✗	✗	✗	✗	✗

On which skin tones you may use the IPL hair removal device:
white, yellow, light brown, brown

On which skin tones you mustn't use the IPL hair removal device:
dark brown, black

On which hair colors you may use the IPL hair removal device:
black, dark brown, brown, light brown

On which hair colors you mustn't use the IPL hair removal device:
light blonde, red, white / gray

1. Use the IPL hair removal device on clean and dry skin.
2. Use a razor to shave the area of the skin on which you want to use the IPL hair removal device. Hair must be shorter than 1 mm. If you use the IPL hair removal device on unshaved skin, it may cause pain, skin redness, or burns. You may also damage the IPL hair removal device this way.
3. Do not use oily substances – such as perfume, lotions, or creams – on the skin that you want to depilate.
4. Make sure that the IPL hair removal device's lamp window is clean. If it isn't, wipe it thoroughly with a dry cloth.

How to use the IPL hair removal device

Use the IPL hair removal device safely

01. Before you use the IPL hair removal device, cover all tattoos and spots on your skin, such as freckles, moles, or birthmarks. Do not use the IPL hair removal device in those places.
02. Do not look directly at the device when it's emitting light.
03. Do not use the IPL hair removal device on nipples or reproductive organs.
04. Do not use the IPL hair removal device on areas of the skin on which you want hair to grow in the future.
05. Do not flash the light from the IPL hair removal device when the device isn't touching your skin. Do it only when the lamp's light exit window touches the skin completely.
06. Do not use alcohol or any other flammable liquids – such as perfume or a disinfectant – to clean the skin on which you'll be using the IPL hair removal device.
07. Using an IPL hair removal device may cause temporary changes in skin color. If you have been using the IPL hair removal device for some time and the color of your skin changed significantly – consult a doctor.
08. Start using the IPL hair removal device in the 1st intensity mode. The higher intensity mode you set, the greater is the likelihood of side effects. If you don't feel any discomfort while using the IPL hair removal device, you may set higher intensity mode next time you use the device.

How to use the IPL hair removal device

01. Remember to shave the area of the skin on which you want to use the IPL hair removal device. You may use the shaving razor from the set.
02. Connect the IPL hair removal device to the power adapter and plug the power adapter in. You'll see the word OFF on the display.
03. Put on the safety glasses which are in the set.
04. Touch and hold the  button for 3 seconds to turn on the IPL hair removal device.
05. The IPL hair removal device will start working in the 1st intensity mode. If you're using the device for the first time, use this mode. If you don't feel any discomfort, you may increase the intensity mode next time you use the device.
06. To set higher intensity mode, touch the  button. The IPL hair removal device has 5 intensity modes. Use the mode that makes you feel slight discomfort while using the IPL hair removal device. Then the device will work most efficiently.
07. Put the IPL hair removal device's lamp window on the skin that you want to depilate. Hold the device at a 90-degree angle to your skin.
08. Press the **M** button to emit light.
09. After you depilate an area of the body, e.g., legs or bikini line, you may repeat it once or twice. Thanks to that, you'll get better results. Do not use the IPL hair removal device more than 3 times in a row on the same area of your skin.
10. You may use the IPL hair removal device in an automatic mode. To turn this mode on, press the **M** button and hold it for 3 seconds. In this mode, the IPL hair removal device emits light automatically. You just need to move it around the area of the skin that you want to depilate. To turn the manual mode on, press the **M** button and hold it for 3 seconds.

Tip | You may use the IPL hair removal device with a cooling function. When this function is turned on, the device cools the skin after it flashes light. To turn this function on, touch the  button.

11. When you've finished using the IPL hair removal device, touch the  button and hold it for 3 seconds to turn the device off.
12. Unplug the power adapter.
13. If you don't turn the IPL hair removal device off, it will turn off automatically after 3 minutes.

Specifications

Product name	IPL Hair Removal Device Jilt
Model	LY1013
Input voltage	100-240V
Output voltage	12 V
Current	4 A
Power	48 W
Display	LCD
Weight	265.3 g

Warranty and technical support

Your product is covered by a 24-month manufacturer's warranty. For more information go to <https://www.lykkele.pl>. If you have questions about using the product – contact us at kontakt@lykkele.pl. Manufacturer: al.to sp. z o.o., ul. Bojemskiego 25, 42-202 Częstochowa, Poland.

Disposal and the EU compliance statement

We, as the manufacturer of this equipment, declare that it meets the rules of the appropriate EU directives. If you need a copy of the EU Declaration of Conformity – contact us.

Do not throw this equipment out with other household waste. This equipment is made from materials that may be harmful to the environment and human health, if the product is disposed of in the wrong way. When you need to throw away an old product, take it to a designated collection point.

Do not use this equipment outside. Do not expose it to various weather conditions.

Do not use this equipment near a bathtub, shower, sink, or any other container with water. If you're using this device in a bathroom, unplug it after you've finished using it. Even a turned-off device may be dangerous when it's near water.

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