LYKKELE


Lien

Irrigator

EN

Lien

Irrigator

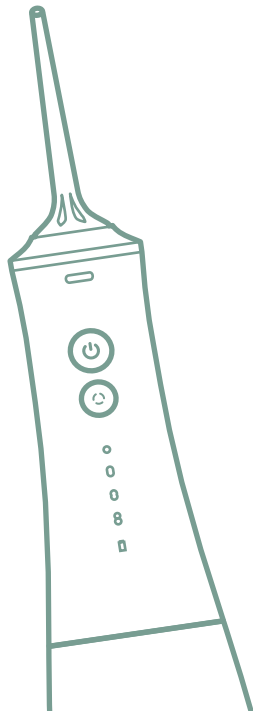Lien

Irrigator

This is Lien, the irrigator from Lykkele. The product will allow you to take even better care of your smile. You will clean your teeth, gums and tongue. See how to use it. **Read its user manual before use. Keep it for future reference.**

Table of contents


01. Safety information
02. Description and casing
03. Irrigator tips
04. Modes of use
05. Who should not use the irrigator
06. Prior to first use of the irrigator
07. How to use the irrigator
08. After use
09. Battery charging and maintenance
10. Most common problems and troubleshooting
11. Specifications
12. Warranty and technical support
13. Environmental protection and compliance with law

01. Safety information


- Use the irrigator in accordance with the user manual.
- Do not use the irrigator if it is damaged - for example, it is not working properly, has a cracked casing, or the plug of the charging cable is damaged.
- Do not modify or repair the irrigator yourself. We shall not be held liable for damage or damage to the product caused by user intervention.
- Keep the irrigator out of reach of children.
- Only use irrigator accessories recommended by the manufacturer.
- If the irrigator falls into the water while charging, disconnect the power plug before removing it from the water.
- Do not press the nozzle against your gums or teeth to avoid their damage.
- Do not use if nozzle is missing.
- Persons suffering from periodontal disease or persons who have undergone oral surgery within the past 2 months should consult a physician before using the product.
- If you drop the irrigator, replace the tip, even if it does not look damaged.
- Unplug the cord before cleaning the irrigator or charging cable. Clean with a slightly damp cloth.
- When unplugging the charger, hold it by the plug, not the cable. Do not touch the plug with wet hands.
- Do not expose the product to fire or high temperatures.
- To avoid burns, do not use water hotter than 40 °C.
- Remember to clean the mouthpiece regularly.
- Some mouthwashes may corrode this product. Do not use them.

02. Description and casing


The irrigator is designed to remove plaque and effectively clean hard-to-reach places in the mouth, such as spaces between the teeth and periodontal pockets. Thanks to a specially generated stream of water under controlled pressure, the irrigator provides thorough cleaning of teeth, gums and tongue. Our irrigator offers a variety of tips and modes that can be adapted to your individual needs. This way, you may choose the most appropriate option to achieve optimal results in your daily oral hygiene.

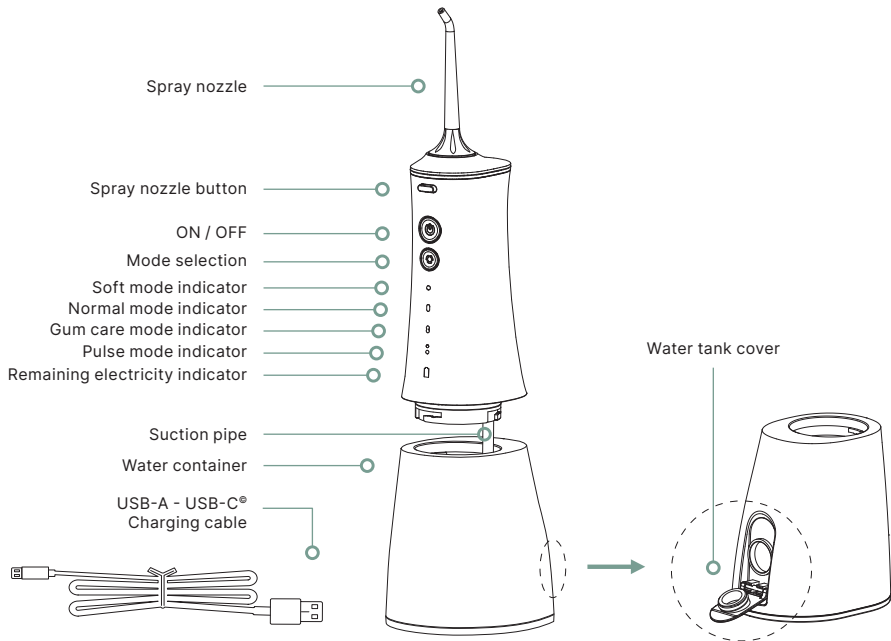
03. Irrigator tips

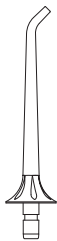A

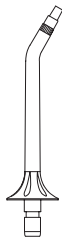B

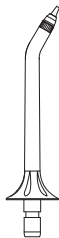C

D

The set accompanying the irrigator includes 6 different tips, each of which performs specific functions and allows effective cleaning of different areas of the oral cavity. Below is the description of individual tips:

- A) Basic head** (2 pieces) - designed for everyday use, especially for healthy teeth and gums. You can effectively remove plaque and food particles from hard-to-reach places, such as interdental spaces and the lower parts of the gums. In addition, by massaging your gums with this tip, you can strengthen their condition.
- B) Tip for braces** (2 pieces) - specially designed to clean areas around the braces. It will remove sediment and food residues from the braces, as well as effectively remove plaque from the enamel. You may also use it to massage your gums.
- C) Tip for cleaning periodontal pockets** - intended for cleaning periodontal pockets and hard-to-reach places in the oral cavity. It can also be used to introduce liquid medications between the gum line and the periodontal pocket. We recommend using it only in a soft mode.
- D) Tongue cleaning tip** - of a small spatula shape, used to remove coating from the surface of the tongue. Use of this tip regularly to maintain proper bacterial flora in the mouth and fresh breath.

04. Modes of use

The irrigator has 4 modes of use:

- A) Soft** - use it if you have sensitive teeth and gums. Use this mode also for the first 2 weeks of use, even if your teeth and gums are not sensitive.
- B) Normal** - use it for daily cleaning.
- C) Gum care** - alternating pulsation of the water stream: strong and weak. For massaging gums and interdental spaces.
- D) Pulse** - pulsating mode for massage. Precise and effective rinsing.

How to change the modes

Press the On/Off button to turn on the irrigator. Select the desired mode by pressing the mode button (D). You can also switch modes while using the product using the mode button (D). The modes change in the order they are listed on the product.

05. Who should not use the irrigator

The irrigator should not be used by:

- Children up to 8 years of age,
- Persons who have ulcers in the mouth,
- Persons who have open sores in the mouth.

Some persons should use the irrigator in the presence of a caregiver. These are:

- Children up to 8 years of age,
- Persons with disabilities if they cannot use the irrigator on their own.

Recommendations

Your oral health is a priority, so we always recommend consulting with a dental or orthodontic professional prior to starting any new oral hygiene practice. This should be followed especially if you are undergoing dental treatment, have dental implants, piercings in the mouth or wear braces.

06. Prior to first use of the irrigator


Recommendations

The gums may bleed during the first two weeks of using the irrigator. To prevent it, we recommend using the irrigator in a soft mode during this period. Regular use should stop the bleeding. If bleeding continues for more than two weeks or is heavy, we recommend consulting a specialist.

Remember!


Fully charge the battery prior to first use of the irrigator. It will take about 5 hours (from 0 to 100%). Do not use the product while it is charging. To charge the irrigator, connect one end of the cord to the product and the other to a 5V 2A wall charger (not included). Plug the charger into an electrical outlet. You can also connect the cable to another power source with appropriate voltage, e.g. power bank or laptop.

01. Insert the tip into the hole at the top of the irrigator and push it until you hear a click

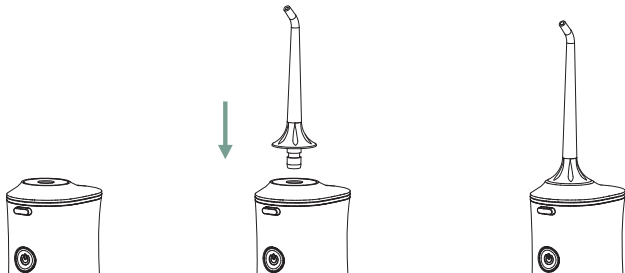
02. Fill in the water container. Use one of the two below methods:

- a) Open the water container cover. Next, holding the product horizontally, pour water into the container and close the cover.

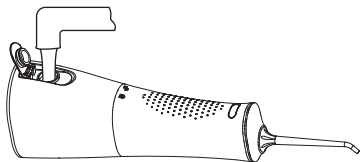
- b) Unscrew the water container to remove it. Hold the water container upright and fill with water. When assembling, press it upwards and turn it towards the lock to lock it.

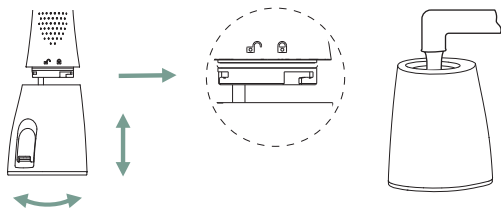Recommendations

Make sure the irrigator is turned off before unscrewing the container and pouring water into it.

03. To remove the tip, press the release button (B).

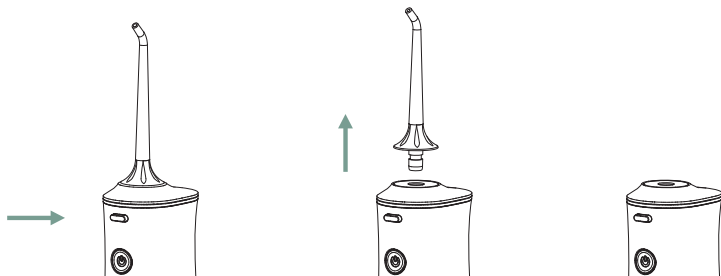07. How to use the irrigator

01. Hold the product diagonally and insert the irrigator tip into your mouth. Point the nozzle at the teeth.



02. Close your mouth to prevent water splashing.
03. Turn on the irrigator by pressing the ON/OFF button.
04. Rotate the irrigator tip to change the direction of the water stream.

Remember

Do not direct the water stream under the tongue.


05. Clean your teeth by moving the irrigator from the back teeth to the front teeth. Direct the water stream along the gum line.
06. To thoroughly clean interdental spaces, make pauses between the teeth.

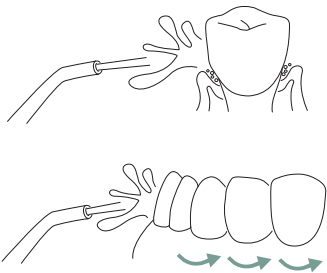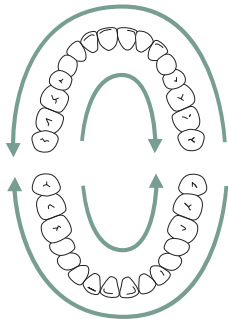07. Clean your teeth with the irrigator for 2 minutes. After this time, the product will automatically turn off.
08. If you want to turn off the irrigator earlier, press the On/Off button.



08. After use

01. Turn off the irrigator if it did not turn off by itself.
02. Open the water container cover.
03. Discard the water from the product.

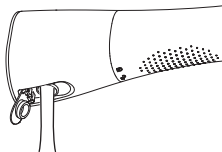09. Charging and battery maintenance

Battery charging

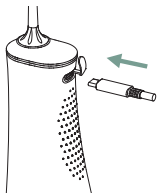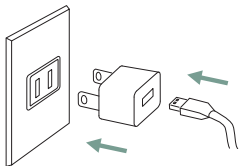01. If you notice a decrease in power when using the irrigator and the battery level indicator starts flashing, it means that the battery is about to run out. Connect the irrigator to the charging device.
02. While the battery is charging, the battery level indicator will flash.
03. When the battery is fully charged, the battery level indicator will be on continuously.

!Backlight!

When the battery level is low, you will see a red flashing light.
While charging, the green LED flashes.
When the device is fully charged, the LED is green.

Battery maintenance

- If you know that the irrigator will not be used for an extended period of time, we recommend you to fully charge the battery beforehand.
- To ensure the longest possible battery life, charge it when you notice a drop in power of the irrigator.

Proper maintenance of the irrigator battery is important to maintain optimal device performance. Ensuring proper charging and responsible management of battery usage will allow you to use the irrigator for a long time.

10. Troubleshooting

Problem	Cause	Solution
Water is leaking from where the irrigator tip is inserted.	The ring that seals the tip insertion point is damaged and the tip cannot be inserted properly.	Contact your retailer to repair the ring.
Water container is leaking.	The seal is loose.	Try to screw in the tank properly once again.
The water pressure is low.	Voltage is too low.	Charge the battery.
Irrigator does not want to turn on.	Voltage is too low.	Charge the battery.

11. Specifications

Product	Irrigator
Model	Lien
Power supply	3.7 V
Power	5 W
Modes	Soft:50-70 PSI • Standard:70-120 PSI • Gumcare:50-80 PSI • Pulse:50-70 PSI
Water pressure	20-120 PSI (1.37 - 8.27 Bar)
Pulses frequency	2000 times per minute
Tips in the set	2 x Basic tip • 2 x Orthodontic braces tip • 1 x Tongue cleaning tip • 1x Tip for cleaning periodontal pockets
Water container capacity	280 ml

12. Warranty and technical support

Your product is covered by a 24-month manufacturer's warranty. For more information go to <https://www.lykkele.pl>. If you have questions about using the product – contact us at kontakt@lykkele.pl. Manufacturer: al.to sp. z o.o., ul. Bojemskiego 25, 42-202 Częstochowa, Poland.

13. Environmental protection and compliance with law

We, as the manufacturer of this equipment, declare that it meets the rules of the appropriate EU directives. If you need a copy of the EU Declaration of Conformity – contact us.

Do not throw this equipment out with other household waste. This equipment is made from materials that may be harmful to the environment and human health, if the product is disposed of in the wrong way. When you need to throw away an old product, take it to a designated collection point.

Do not use this equipment outside. Do not expose it to various weather conditions.

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