
**LYKKELE**

**Berg**


Massager

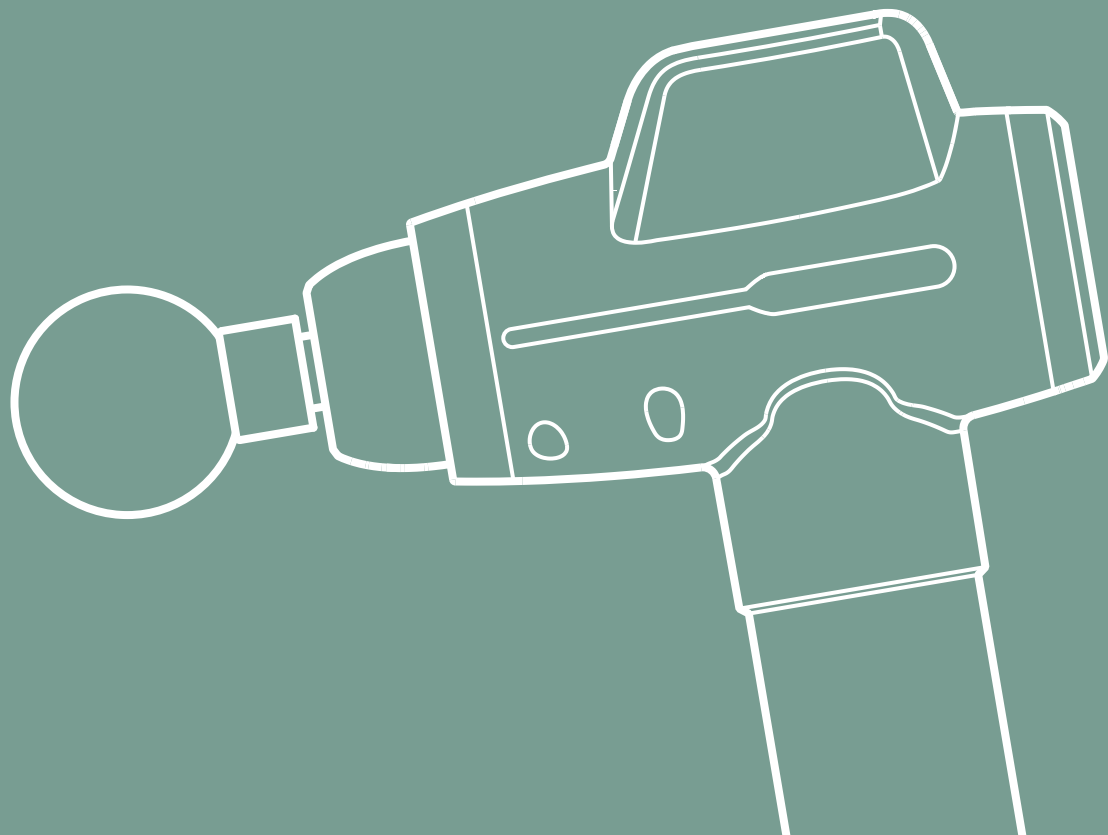EN

# Berg

## Massager
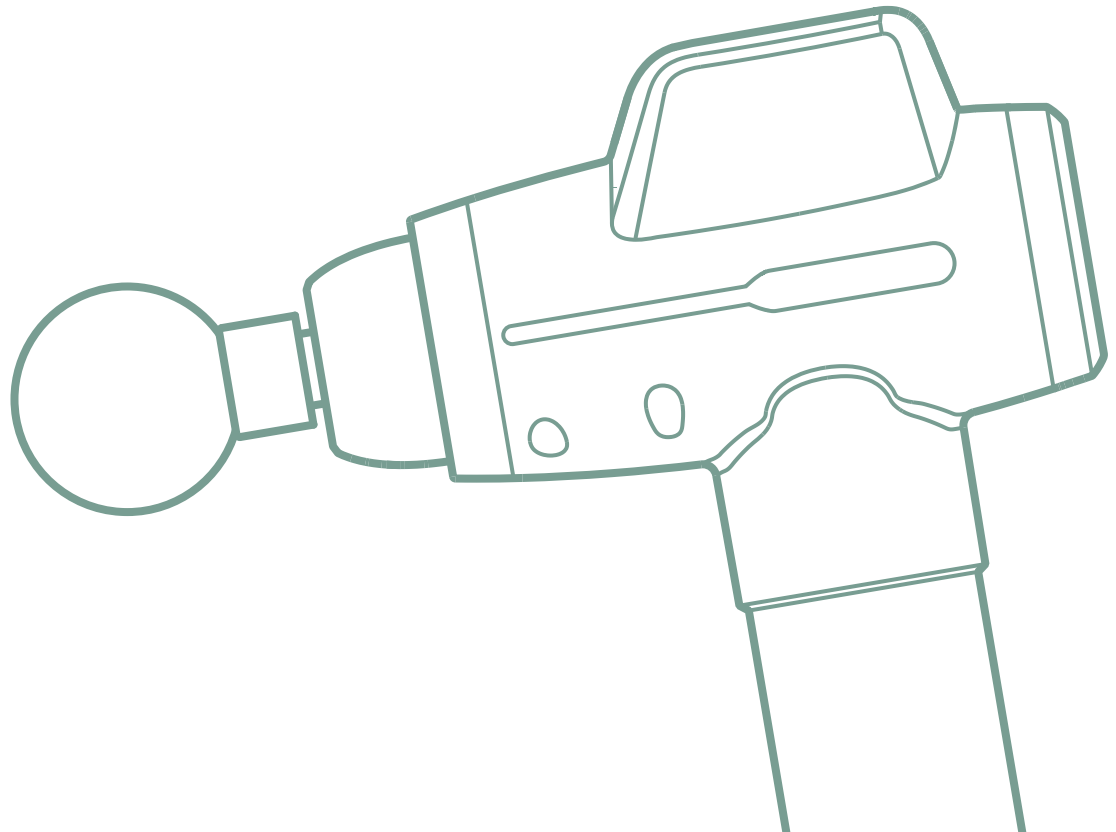
# Berg

## Massager

Here's the Berg - a massage gun designed to relieve muscle tension and bring your body back into balance. See how to use it.  
**Before using the product, read the instruction manual. Keep it for future reference.**

## Table of contents

01. Safety information
02. Description and structure
03. Heads
04. How to use the massage gun
05. Specifications
06. Warranty and technical support
07. Disposal and the EU compliance statement

## 01. Safety information

Persons who meet any of the following conditions should consult a doctor before using a massage gun:

- Pregnant women
- Persons with pacemaker
- Persons who have recently undergone surgery
- Persons with spinal, joint or bone defects
- Persons with complications caused by other diseases (such as neuropathy or retinal damage)


**If you are unsure whether to use a massage gun - consult your doctor or specialist.**

- The product is intended for healthy adults only. Children are not allowed to use a massage gun.
- Use the product only on soft tissues, avoiding causing pain or discomfort.
- Do not use on areas of the body where there are fewer muscles, e.g.: head, tibia, back of feet, back of palms, etc.
- Regardless of how you use it, bruising on the body may occur. Check the massage area regularly and if pain or discomfort occurs, stop using the product immediately.
- Use only on a dry and clean skin surface, gently pressing and gliding over the skin.
- Keep fingers, hair and other body parts away from the shaft and vibrating head.
- Do not use around sensitive body parts such as eyes, teeth, enlarged breasts, reproductive organs, prostheses, etc.
- Do not use the product on areas where you suspect a fracture has occurred.
- If the pain is persistent, stop using the massage gun immediately and consult your doctor. This product is not a medical device and is not intended for medical purposes.
- Use a suitable vibration head for optimum effect.
- Do not block or place any objects in the ventilation openings of the product.
- Avoid pouring water and getting water in the ventilation openings.
- Check the condition of the product carefully before use. If in doubt, do not use it and contact us.
- Do not make any modifications to the product. Do not repair the unit yourself.
- Do not leave the massage gun switched on or connected to the charger unattended.
- Do not start the device while it is charging.

- Do not throw the device into fire.
- Do not charge batteries near a fire or in extreme heat conditions. Avoid use near heat sources such as fire or radiator.
- If the battery leaks or emits an unpleasant odour, immediately place the unit away from a heat source.
- Do not throw the device into fire or heat it up.
- Do not store together with or transport in the presence of other metal objects.
- Do not pierce the battery case with any sharp objects, do not throw or hit.
- Do not place the unit in a microwave oven.
- Fully charge the battery before first use (approximately 3-4 hours).
- You can charge the battery at any time, however, make sure that the massage gun is switched off during charging.
- We do not recommend discharging the battery completely.
- The average working time of the massage gun is approximately four hours, but this is influenced by the frequency of vibration and the pressure applied.


## 02. Description and structure

The massage gun is designed to relieve muscle tension and effectively relax hard-to-reach areas on the body, such as the muscles of the back or legs. Thanks to specially generated vibrations, it ensures thorough relaxation of muscles, connective tissue and skin. Our massage gun offers different heads and modes that can be customised to suit your personalised needs, ensuring optimal results in your daily body care routine. Remember to use the massager as instructed, avoiding areas of the body that may be sensitive.

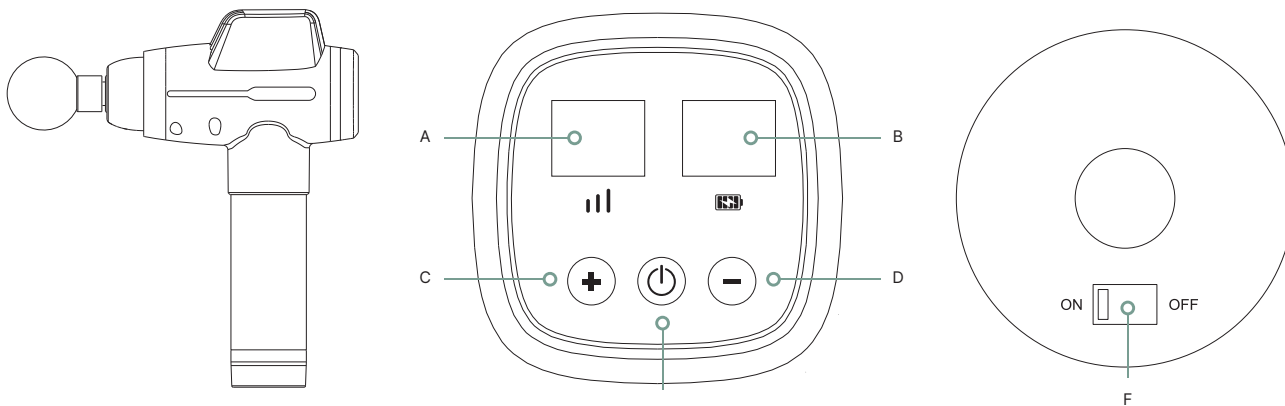A	Frequency Displaying: Shows the current vibration frequency	D	Speed - : Push button to decrease vibration speed
B	Power Display: Indicates battery charge level	E	Start/Stop: Push button to start and stop the massage
C	Speed + : Push button to increase vibration speed	F	On/Off: switch on the bottom of the unit. When it is in

## 03. Heads

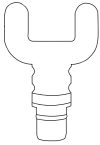
### U-shaped head:

perfect for massaging the muscles around the spine, such as the neck and lumbar. Its unique shape allows it to reach sensitive areas safely and effectively.


### Cylindrical head:

designed for deep point massage. Its narrow, cylindrical shape makes it possible to reach deep muscles and release points.

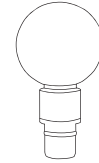
### Ball head:

designed to massage large muscle groups such as the back or legs. Its circular shape allows for wide coverage and effective muscle relaxation.
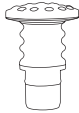
### Flat head:

designed for general massage. Its shape allows for wide coverage and relief of tension in different parts of the body.
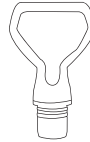
### Spade head:

designed to massage the areas around joints and bones. Its shape allows it to reach these areas, protecting the sensitive bone structures.

### Mushroom head:

designed to massage large muscle groups such as back or thighs. Its wide, rounded surface allows for gentle and sensitive massage of muscles.

Tip

Switch off the unit to install the massage head. Insert the head with the appropriate side into the opening, then press it down. To remove the head, pull it in the opposite direction.

## 04. How to use the massage gun

01. Select an appropriate head and install it .
02. Switch on the unit by sliding the switch (F) on the bottom of the handle.
03. Press the Start/Stop button (E) to start the massage gun.
04. Adjust the frequency level using the Speed + (C) / Speed - (D) keys from 1 to 30.
05. If "E1" message appears in the battery window, it means that the battery is discharged and needs to be charged.

Charging levels		
E1 – 0%	1 - 10%	10 - 100%

06. The massage gun is equipped with a timer function and automatically switches off after 15 minutes of operation.
07. If you press the head too hard, the massager will trigger the protection function and automatically switch off (message E3). After switching off the massager, you can restart it immediately.
08. Once you have finished using the device, switch it off using the push button(6) on the bottom of the handle. If the massager is not used for 3 minutes, it will automatically enter standby mode. It can be restarted by pressing any button or moving the switch on the bottom of the handle.

### Example of a massage technique:

Make sure the gun is fully charged before starting the massage. Select the appropriate head depending on the area you wish to massage and attach it to the gun.

Area	Head	Duration
Spinal muscles	U-shaped head	2 minutes
Back muscles	Ball head	3 minutes
Abdominal muscles	Mushroom head	2 minutes
Arm muscles	Flat head	2 minutes (1 minute for each arm)
Leg muscles	Ball head	4 minutes (2 minutes for each leg)
Muscles around the joints	Shovel-shaped head	2 minutes

### Massage technique:

Guide the massage gun along the muscles, using movements from the inside to the outside. Press the gun lightly against the body, but do not use too much pressure to avoid injuring the muscles. Remember to breathe deeply during the massage to help the muscles relax.

## 05. Specification

Product	Massage gun
Model	Berg
Amplitude	12 mm
Rated power	22.2 W
Rated voltage:	12 V
Battery capacity	2500 mAh
Charging time	3-4 hours
Working time	Approximately 4 hours (depends on vibration frequency and pressure applied)
Accessories included	• Power supply • U-shaped head • Cylindrical head • Ball head • Flat head • Shovel-shaped head • Mushroom head
Weight	740 g (without head)=
Dimensions	172 mm (without head) x 257 mm x 67.5 mm

## 06. Warranty and technical support

Your product is covered by a 24-month manufacturer's warranty. For more information go to <https://www.lykkele.pl>. If you have questions about using the product – contact us at [kontakt@lykkele.pl](mailto:kontakt@lykkele.pl).  
 Manufacturer: al.to sp. z o.o., ul. Bojemskiego 25, 42-202 Częstochowa, Poland.

## 07. Disposal and the EU compliance statement

We, as the manufacturer of this equipment, declare that it meets the rules of the appropriate EU directives. If you need a copy of the EU Declaration of Conformity – contact us.

Do not throw this equipment out with other household waste. This equipment is made from materials that may be harmful to the environment and human health, if the product is disposed of in the wrong way. When you need to throw away an old product, take it to a designated collection point.

Do not use this equipment outside.  
 Do not expose it to various weather conditions.

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