# LYKKELE

## Edda

Sonic Toothbrush

EN

# Edda

Sonic Toothbrush

# Edda

## Sonic Toothbrush

This is Edda, Lykkele sonic toothbrush. Thanks to its powerful vibrations, it will clean your teeth more precisely than a manual toothbrush. Your smile will be even whiter with its help. Check out how to use it. **Before you use the sonic toothbrush, read the manual carefully. Do not throw the manual out.**

## Table of contents


01. Safety information
02. The sonic toothbrush's description and parts
03. The sonic toothbrush's tips
04. The sonic toothbrush's modes
05. Who shouldn't use the sonic toothbrush
06. Before you use the sonic toothbrush for the first time
07. How to use the sonic toothbrush
08. How to clean the sonic toothbrush and the charger
09. Charging and battery maintenance
10. Troubleshooting
11. Specifications
12. Warranty and technical support
13. Disposal and the EU compliance statement

# 01. Safety information

- Use the sonic toothbrush only as we indicate in this manual.
- Do not turn the sonic toothbrush on if it's damaged – e.g., it doesn't work properly, its case or its charger's plug is broken.
- If the sonic toothbrush is damaged, send it to us for repair. Otherwise, your warranty expires. We are not responsible for any damage that a sonic toothbrush which wasn't repaired in our service may cause.
- Keep the sonic toothbrush out of children's reach.
- Only use tips and accessories that the manufacturer recommends.
- Do not use the sonic toothbrush nor the charger outdoors.
- You may use the sonic toothbrush in temperatures from -5 up to 50 °C.
- Do not expose the sonic toothbrush nor the charger directly to a heat source or sunlight.
- If you're not using the charger - unplug it.
- Before you clean the sonic toothbrush or its charger, unplug the charger. Clean the device and the charger with a slightly damp cloth.
- When you unplug the charger, hold onto the plug, not the cord. Do not touch the plug with wet hands.

# 02. The sonic toothbrush's description and parts

The sonic toothbrush is a device that will help you thoroughly clean your teeth, gums and tongue. If you regularly use this toothbrush to clean your teeth, it will help you remove dental plaque. It will also help you prevent the accumulation of dental plaque on your teeth. With this toothbrush, you'll be able to remove teeth stains caused by drinking coffee or tea. The device will help you maintain a healthy oral cavity. You may use different modes to adjust the toothbrush to your needs.
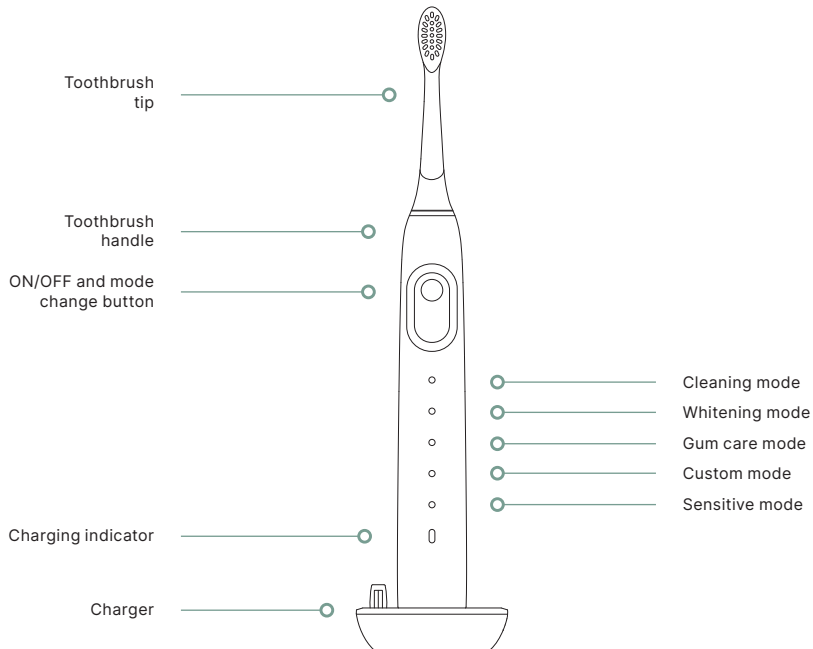
## 03. The sonic toothbrush's tips

You'll find 2 basic toothbrush tips in the set. Replace the toothbrush tip every 2-3 months to clean your teeth best and to maintain healthy oral cavity.

## 04. The sonic toothbrush's modes

### The toothbrush has 5 modes:

- **cleaning** - use it to brush your teeth daily.
- **whitening** - use it to remove teeth discoloration.
- **gum care** - use it to massage your gums. If you do it regularly, you'll improve their blood circulation. It will help you maintain healthy gums.
- **custom** - in this mode, you may set the intensity of toothbrush vibrations. Turn this mode on, then press and hold the ON/OFF button. Toothbrush vibrations will change – you may choose out of 10 intensity levels. When the vibrations suit you, release the button. The toothbrush will record this setting.
- **sensitive** - choose this mode if you've never used a sonic toothbrush before. Use it also if you have sensitive teeth or gums.

### How to change modes

Press the ON/OFF button to turn the toothbrush on. When you turn it on for the first time, it starts working in the **cleaning** mode. To enter the function that changes modes, press the ON/OFF button twice. Once the function that changes modes is on, press the ON/OFF button to switch between modes. Modes change in the order in which they are listed on the device.

### The toothbrush remembers the mode

If you use a mode for more than 5 seconds, the toothbrush will remember it. Next time you turn the toothbrush on, it will work in this mode.

## 05. Who shouldn't use the sonic toothbrush

Who shouldn't use the sonic toothbrush:

- kids younger than 12;
- people who have loose dentures;
- people who have bleeding gums; and
- people who have hypersensitive teeth or gums.

**Tip** | People with disabilities should use the sonic toothbrush with a guardian present if they can't use it on their own.

**Tip** | Consult your dentist whether you may use the sonic toothbrush if:

- you've had periodontal surgery less than two months ago; or
- your gums are bleeding after you brush your teeth with the sonic toothbrush, and it lasts more than a week.

**Tip** | If you have a pacemaker or other implanted electronic device, consult your doctor if you may use the sonic toothbrush.

## 06. Before you use the sonic toothbrush for the first time

**Tip** | When you use the sonic toothbrush for the first time, you may feel irritation or numbness in the oral cavity. These symptoms are normal. They will gradually disappear when you use the sonic toothbrush for longer than two days. The symptoms are usually completely gone after a week.

### Charge the sonic toothbrush's battery

Before you use the toothbrush for the first time, fully charge its battery. It will take about 18 hours.

To charge the toothbrush's battery:

01. Plug the charger's cord into a 5V 2A USB charger.
02. Then plug the USB charger in the electric socket. You may also plug the charger's cord into another 5V 2A power source, e.g., a power bank or a laptop.
03. Place a turned-off toothbrush on the protrusion in the middle of the charger.
04. When the toothbrush is charging, do not press the ON/OFF button.

### Insert the toothbrush tip on the toothbrush handle

01. Insert the toothbrush tip on the metal shaft on top of the toothbrush handle. Align the toothbrush tip so the bristles point the same direction as the front of the toothbrush handle.
02. Press the toothbrush tip down onto the metal shaft until it stops.

## 07. How to use the sonic toothbrush

01. Wet the toothbrush bristles.
02. Apply toothpaste on the bristles.
03. Place the toothbrush tip inside your mouth and put it against your teeth at a 45 degrees angle – aim it towards your gums.
04. Press the ON/OFF button to turn the toothbrush on.
05. Choose a mode. When you turn it on for the first time, it starts working in the **cleaning** mode.


**Tip** | If it's your first time using a sonic toothbrush or you have sensitive teeth or gums, start with the **sensitive** mode. To turn this mode on, press the ON/OFF button twice – you'll enter the function that changes modes. Then press the button four times to turn on the **sensitive** mode.

06. Start cleaning from the back teeth and move toward the front ones. It will be best if you divide your teeth into 4 areas and then clean these areas one after another.

10 seconds for outside of the teeth

10 seconds for topside of the teeth

10 seconds for the inside of the teeth

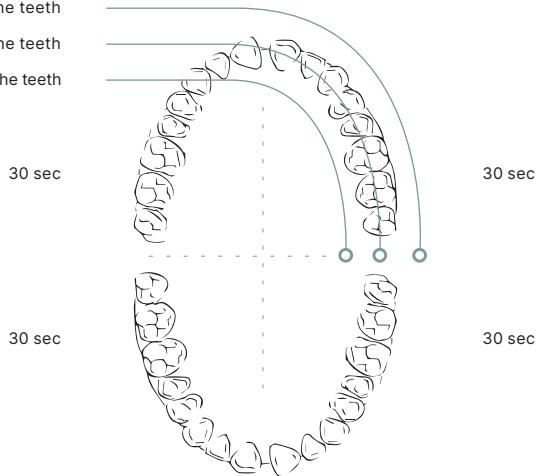07. The toothbrush will pause every 30 seconds to remind you that you should change the area which you're cleaning.
08. Brush your teeth with the sonic toothbrush for 2 minutes. After this time, the device will turn off automatically.
09. If you want to turn the toothbrush off earlier, press the ON/OFF button. The toothbrush must be working in the same mode for at least 5 seconds to turn it off.

## 08. How to clean the sonic toothbrush and the charger

### How to clean the toothbrush tip

01. Rinse the toothbrush tip with water after each use.
02. Do not pull the bristles out of the toothbrush tip.
03. Replace the toothbrush tip every 2-3 months.

### How to clean the toothbrush handle

01. Clean the toothbrush when it's turned off.
02. Remove the toothbrush tip. To do it, pull it up.
03. Rinse the metal shaft with water and then dry it with a cloth. Do it once a week.
04. Clean the toothbrush handle with a soft damp cloth.
05. Do not put the toothbrush under water for longer than half an hour.

### How to clean the charger

01. Each time that you want to clean the charger, unplug it first.
02. Clean the charger with a soft damp cloth.
03. Clean the charger's plug with a dry cloth regularly. If the plug is covered in dust, the charger won't work properly.
04. Do not put the charger under water or any other liquid.

## 09. Charging and battery maintenance

### Battery charging

- When you fully charge the battery before you use the toothbrush for the first time, keep using the toothbrush until the battery dies. Then fully charge it. After this, you may charge it when you want.
- When there's less than 10% of the battery left, the indicator glows red – you should charge the toothbrush.
- When the battery is charging, the indicator flashes green. It takes about 18 hours to fully charge the battery.
- When the indicator glows green, the battery is fully charged.
- If you place a turned-on toothbrush on the charger, the device will automatically turn off and start charging.
- The toothbrush may heat up when it's charging – it's normal. Its temperature shouldn't exceed 40 °C.

### Battery maintenance

- If you know that you won't be using the sonic toothbrush for a long time, fully charge its battery beforehand. Then charge it once every three months.
- In order to ensure the battery's long life, charge it when the indicator glows red.

## 10. Troubleshooting

Problem	Cause	Solution
Vibrations are weak.	You press the toothbrush on the teeth too hard.	Try to brush your teeth softer.
The toothbrush works shorter than before.	The battery is low.	Make sure that you place the toothbrush on the charger correctly. Charge the toothbrush longer.
The toothbrush doesn't turn on.	The battery is low.	Charge the battery. If that doesn't help, contact us.

## 11. Specifications

Product	Sonic Toothbrush
Model	Edda
Voltage	3.7 V
Power	2 W
Modes	• Cleaning • Whitening • Gum Care • Custom • Sensitive
Vibration frequency	42000 strokes per minute
Tips included in the set	Basic tip × 2
Battery capacity	2000 mAh
Weight	174.9 g

## 12. Warranty and technical support

Your product is covered by a 24-month manufacturer's warranty. For more information go to <https://www.lykkele.pl>. If you have questions about using the product – contact us at [kontakt@lykkele.pl](mailto:kontakt@lykkele.pl). Manufacturer: al.to sp. z o.o., ul. Bojemskiego 25, 42-202 Częstochowa, Poland.

## 13. Disposal and the EU compliance statement

We, as the manufacturer of this equipment, declare that it meets the rules of the appropriate EU directives. If you need a copy of the EU Declaration of Conformity – contact us.

Do not throw this equipment out with other household waste. This equipment is made from materials that may be harmful to the environment and human health, if the product is disposed of in the wrong way. When you need to throw away an old product, take it to a designated collection point.

Do not use this equipment outside. Do not expose it to various weather conditions.

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